

## Ready to Use

## ChromoGel™ Strep B Agar (GBS)

The medium use for isolation and identification of Streptococcus B group bacteria (Opaque background)

Formula	g/l
Peptone, Yeast extract	20
Growth factor	16
Salts	8.5
Chromogenic mix	0.5
Agar	15
Selective mix	8
pH: 7.3 ± 0.2 at 25°C	

## Preparation

The medium has been arranged in sterile petri dish with 22-24ml of medium for ready to use.

Primary packaging: cellophane wrapping

Secondary packaging: cardboard box

## Quality Control acc. to ISO 11133

- Productivity control at 37 ± 2°C / 18-24H

Organism	Inoculated CFU	Result
S.agalactiae ATCC 13813	10 <sup>3</sup> - 10 <sup>5</sup>	Purple colonies >30 CFU
Enterococcus faecalis ATCC 29212	10 <sup>3</sup> - 10 <sup>5</sup>	Dark blue colonies >30 CFU

- Selective control at 37 ± 2°C / 18-24H

Organism	Inoculated CFU	Result
Candida albicans ATCC 10231	10 <sup>6</sup> - 10 <sup>8</sup>	Inhibited
E. coli ATCC 25922	10 <sup>6</sup> - 10 <sup>8</sup>	Inhibited

## Storage

Store the medium tightly closed in a dry place at room temperature and use before the expiry date on the label.

Store temperature: 2-20°C

## Ordering Infomation

Code	Products	Package
CGGBS500	ChromoGel™ Strep B Agar Base	500g/Bottle
CGGBS10S	ChromoGel™ Strep B Supplements	10vials/Box
CGGBS9010	ChromoGel™ Strep B Agar (90mm)	10Plates/Box

## Ready to Use

## ChromoGel™ Strep B Agar (GBS)

Môi trường dùng nuôi cấy và định danh Streptococcus nhóm B (Nền thạch trắng sữa)

Thành phần	g/l
Peptone, Yeast extract	20
Growth factor	16
Salts	8.5
Chromogenic mix	0.5
Agar	15
Selective mix	8
pH: 7.3 ± 0.2 at 25°C	

## Chuẩn bị

Môi trường đã được chuẩn bị trong đĩa petri tiệt trùng với 22-24ml môi trường.

Bao bì đóng gói chính: màng cellophane chuyên dụng

Bao bì phụ : hộp giấy các tông.

## Đánh giá chất lượng theo ISO 11133

- Đánh giá hiệu năng ở 37 ± 2°C / 18-24H

Chủng vi nấm	Nồng độ cấy CFU	Kết quả
S.agalactiae ATCC 13813	10 <sup>3</sup> - 10 <sup>5</sup>	Màu tím hoa cà >30 CFU
Enterococcus faecalis ATCC 29212	10 <sup>3</sup> - 10 <sup>5</sup>	Màu xanh đậm >30 CFU

- Đánh giá chọn lọc ở 37 ± 2°C / 18-24H

Chủng vi nấm	Nồng độ cấy CFU	Kết quả
Candida albicans ATCC 10231	10 <sup>6</sup> - 10 <sup>8</sup>	Không mọc
E. coli ATCC 25922	10 <sup>6</sup> - 10 <sup>8</sup>	Không mọc

## Bảo quản:


Bảo quản môi trường ở nơi khô ráo, tránh ánh sáng và sử dụng trước hạn dùng.

Nhiệt độ bảo quản: 2-20°C


## Thông tin đặt hàng

Code	Products	Package
CGGBS500	ChromoGel™ Strep B Agar Base	500g/Bottle
CGGBS10S	ChromoGel™ Strep B Supplements	10vials/Box
CGGBS9010	ChromoGel™ Strep B Agar (90mm)	10Plates/Box

## ChromoGel™ Strep B Agar

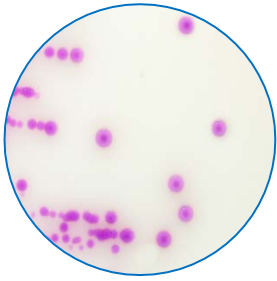

Streptococcus agalactiae

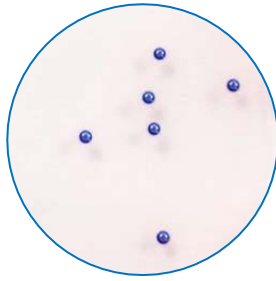

Enterococcus faecalis